



Voyager III <sup>NEW</sup>

 **Belmont**

The Beauty of Dentistry

# FREE 5 Year Warranty

for Voyager Dental Chairs and Dental Units, and Dental Operating Lights.  
Visit [belmontdental.co.uk](http://belmontdental.co.uk) or refer to the back cover for further details.

## Voyager III

### Key improvements

Improved vacuum air-flow  
Chair base and delivery arm re-designed  
Chair ergonomics and comfort improved  
700 Series light (720) with touchless sensor available as standard

### Key features

Ambidextrous  
Ideal console arrangement for Assistant's 'Preparation zone'  
Below-the-Patient swing-arm delivery system  
Mobile delivery systems options  
Seamless upholstery  
Standard or Ultrasoft upholstery colour options  
Air operating system  
Upgradeable Operator's Console functionality

### Unit-mounted Operating Light options

720SLR or 300 LED

### Upgrading and accessorising

Belmont's range of factory built-in options can be used to increase the functionality of a treatment centre and the services that your Practice offers to the Patients. Visit [belmontdental.co.uk](http://belmontdental.co.uk) or see our price list for further details.

### Specification

135kg max. Patient weight | 480mm min. seat height | 780mm max. seat height


Very versatile and compact ambidextrous treatment centres,  
designed and built for the most demanding working environments.



Voyager III **ambidextrous treatment centres**

# Understated versatility & value


Voyager III's versatile Below-the-Patient delivery system and ultra-slim backrest design eliminate physical stresses associated with some treatment centres.

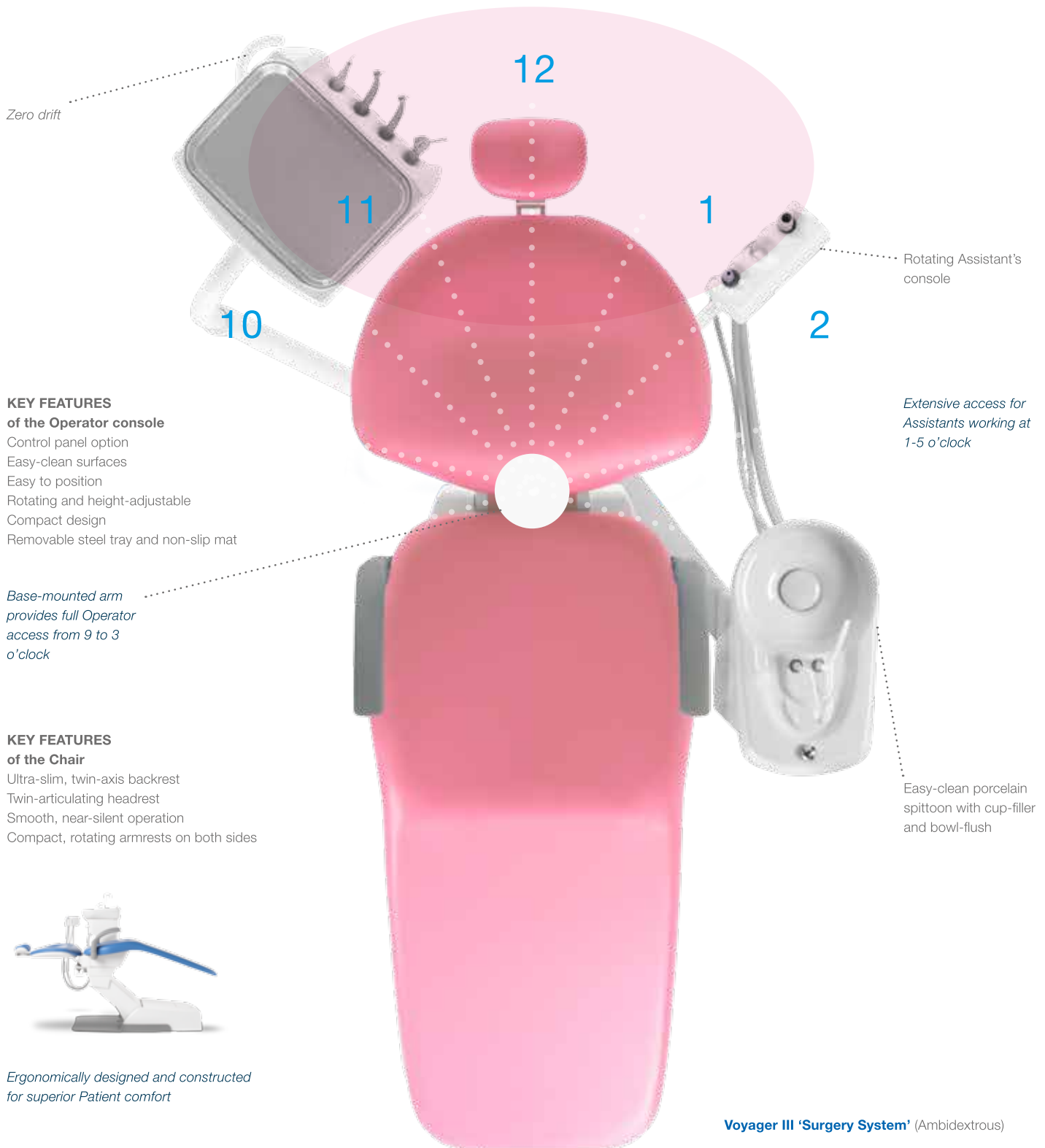


# Left- or right- anytime

## Switchover in no time at all


The central pivot design enables the very easy changeover procedure and provides 100% equality for left- and right-handed Operators.

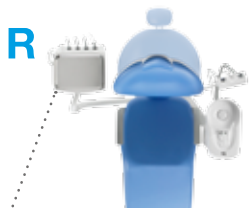


**Right-handed set up.**  
To begin the switchover, raise chair height.

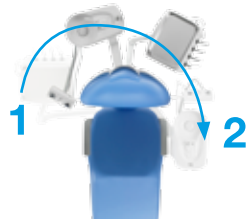

Incline chair.  
**1-2. Swing table round to left side.**

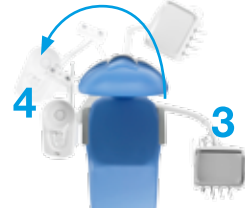

**3-4. Swing cuspidor round to right side.**

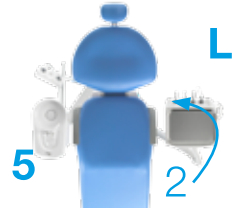

**5. Secure cuspidor in position.**  
**Secure cuspidor locking button.**  
2. Rotate console.  
**Left-handed set up.**



Be more welcoming to every Patient

Obstruction-free zone for easy entry and exit

*Compact, rotating armrests and a low chair position ensure a highly accessible entry*

# Below-the-Patient advantages

Ready to welcome the next Patient, instruments are out of sight and preparation has begun.



Operator's console is discreetly positioned behind the chair in the Assistant's 'Preparation Zone'

## Preparation zone

*Cleaning and preparation can proceed during consultation >*



Easy-clean vacuum filter with detachable tubings

Hygiene-maintenance is easy with the consoles and light gathering together compactly behind the chair.



# Be more efficient

Intelligently designed for a comfortable and compact, easy-clean working environment.



Compact, versatile delivery system offers useful working positions all the way around



## Consultation zone

*An efficient workflow can proceed calmly during cleaning and preparation processes >*

Our advancement of Below-the-Patient delivery system ergonomics realises a greater efficiency of time, space and energy. A level of console versatility yet unseen. Your every task becomes easier; less reaching, a new calmness and clear lines of communication. This is what we strived for, a better daily operation for all.

### Compactness

*Voyager III's footprint is remarkably compact, sparing valuable space, demanding less energy to operate and improving communication at all stages of surgery*

# Control is simple

Control chair movements and instruments at the foot-controls. Voyager III is an uncluttered, simple and reliable working partner which allows the Operator full control with the option of mounting a control panel on the Operator's console.



<sup>1</sup> Requires multi-function control panel (optional)



# Flexible options to suit

**Voyager III 'SP'** ambidextrous chair and cuspidor



Operator's control panel offers two presets, Last-Position Memory, Auto-Return and manual control for chair movements (optional)



**Voyager 'MC'**  
Mobile Cart



**Voyager 'M'**  
Module

# 720SLR for Voyager III

## KEY FEATURES

- Unit-mounted for Voyager III
- White natural beam (4,200 Kelvin)
- 3 light settings:
  - Low (18,000 Lux), High (28,000 Lux) and Composite Curing Safe Mode (at 8,000 Lux)
- 650mm focal distance
- Touchless sensor switch
- Patient mirror
- Auto ON/OFF with preset chair positions
- Halogen bulb
- Twin-axis head movement



**650mm focal distance** and the very compact light head allow extensive working space for the dental team

Light intensity button and indicator: soft-push-up button is easily and hygienically activated with the back of a hand



# 320 LED for Voyager III

## KEY FEATURES

- Unit-mounted for Voyager III
- White natural beam (5000 Kelvin)
- Stepless intensity control:  
3,100 to 28,000 Lux
- 650mm focal distance
- Composite Curing Safe Mode (at 4,300 Lux)
- Touchless sensor switch
- Auto ON/OFF with preset chair positions\*
- 5 LEDs
- Triple-axis head movement



Unique stepless control allows adjustment to any light intensity between 3,100 and 28,000 Lux



# **5**<sup>years</sup> Our warranties are extended and **FREE**

On condition that the purchaser returns the Extended Warranty form to Takara Belmont (UK) Ltd, the Standard Warranty will be extended to ~ **5 Years** for Voyager Dental Chairs and Dental Units, and Dental Operating Lights. **2 Years** for Intra-Oral X-Ray equipment.

**Standard Warranty Policy** – All Belmont products are guaranteed against mechanical, electrical and workmanship faults, under normal working conditions, for a period of one year from date of delivery. Paintwork, plating, natural wear/tear of all tubings are not warranted unless they are proved to be defective at time of delivery. OEM factory built-in options (i.e. scalers, motors and CVS kits) carry a warranty as provided by the individual manufacturer.



[belmontdental.co.uk](http://belmontdental.co.uk)

@BelmontDentalUK

London showroom 020 7515 0333

Manchester showroom 0161 743 9992

All showroom visits by appointment only

The manufacturers reserve the right to alter the specifications without prior notice.

Ref\_BelmontUKWireforks\_2015.10.

